


Delivering Optiflow™ Nasal High Flow

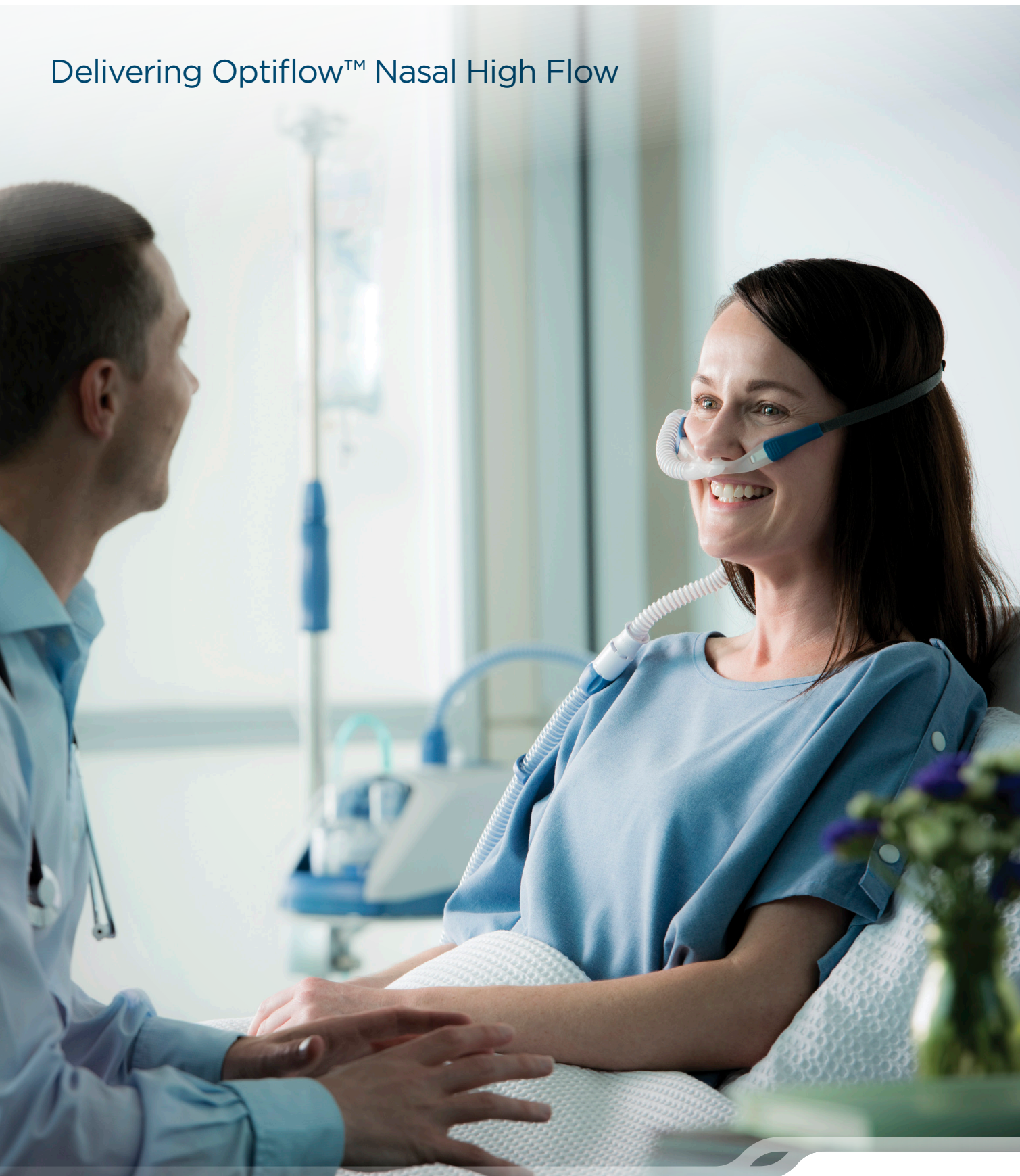

Introducing the F&P AIRVO™ 2, delivering respiratory support and airway hydration with a focus on comfort and ease of use.

The F&P AIRVO 2 has been designed specifically to be an integrated solution for delivering Optiflow™ Nasal High Flow. Fisher & Paykel Healthcare's leading humidification technology allows the AIRVO 2 to comfortably provide high flows of air/oxygen mixtures to spontaneously breathing patients, through the unique Optiflow nasal interface.

Adjustable temperature and flow settings

- Three temperature settings (37, 34, 31 °C) help achieve comfort and compliance
- Integrated flow generator delivers a wide flow range 10 - 60 L/min for adult patients and 2 - 25 L/min for pediatric patients. No wall air supply required

Optiflow™+ Nasal Interface

- Soft, flexible prongs
- Wide-bore design reduces gas jetting
- Prongs contoured to the patient's nose



Supplementary oxygen (when required)

- Oxygen can be added from a wall supply or cylinder
- Inbuilt ultrasonic oxygen analyzer requires no calibration, service or replacement

Airspiral™ heated breathing tube


- Dual spiral heater wires and unique integrated temperature sensor
- No separate temperature probes or heater-wire adapters required

Designed for simple setup, use and cleaning

- Helpful onscreen animations assist with setup and troubleshooting

An integrated system for delivering Optiflow Nasal High Flow

- The F&P AIRVO 2, with its wide range of accessories, provides versatility, mobility and convenience

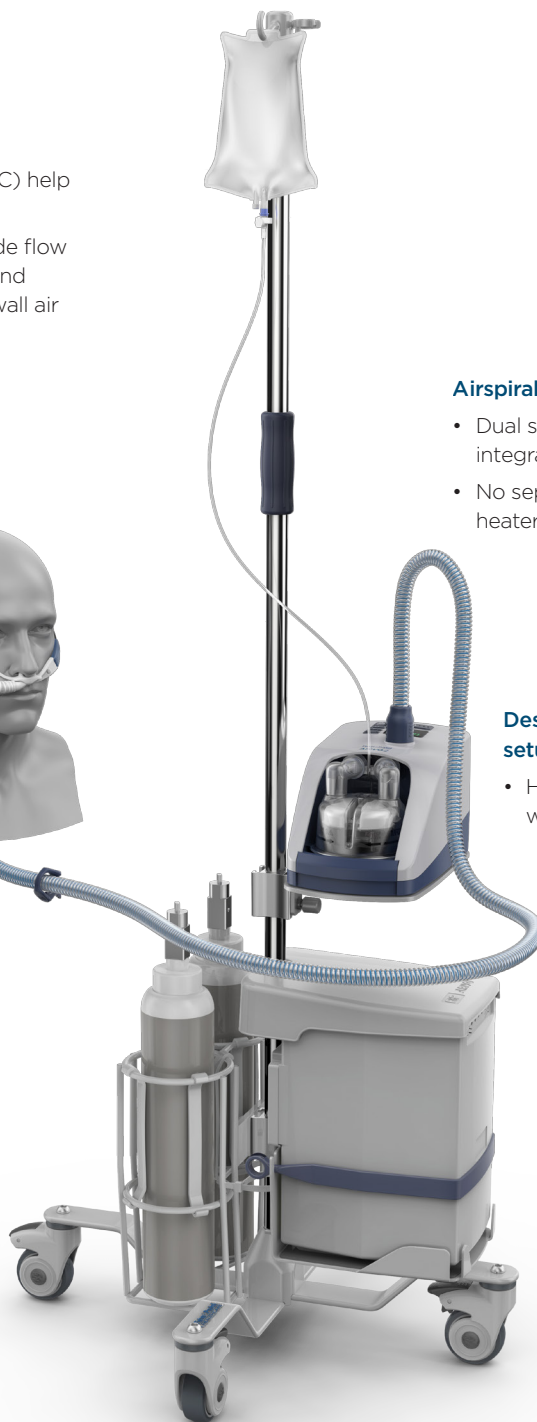

Product codes	
PT101xx	AIRVO 2 humidifier
OPT942	Optiflow+ small interface (20-pack)
OPT944	Optiflow+ medium interface (20-pack)
OPT946	Optiflow+ large interface (20-pack)
OPT316	Optiflow Junior interface - infant (20-pack)
OPT318	Optiflow Junior interface - pediatric (20-pack)
OPT970	Optiflow+ Tracheostomy interface (20-pack)
OPT980	Optiflow+ mask interface adapter (20-pack)
900PT561	AirSpiral heated breathing tube & chamber kit (10-pack)
900PT405	Pole mounting tray
900PT421	Mobile pole stand
900PT426	Plastic basket
900PT427	Oxygen bottle holder
900PT411	UPS mounting kit